


Miscellaneous

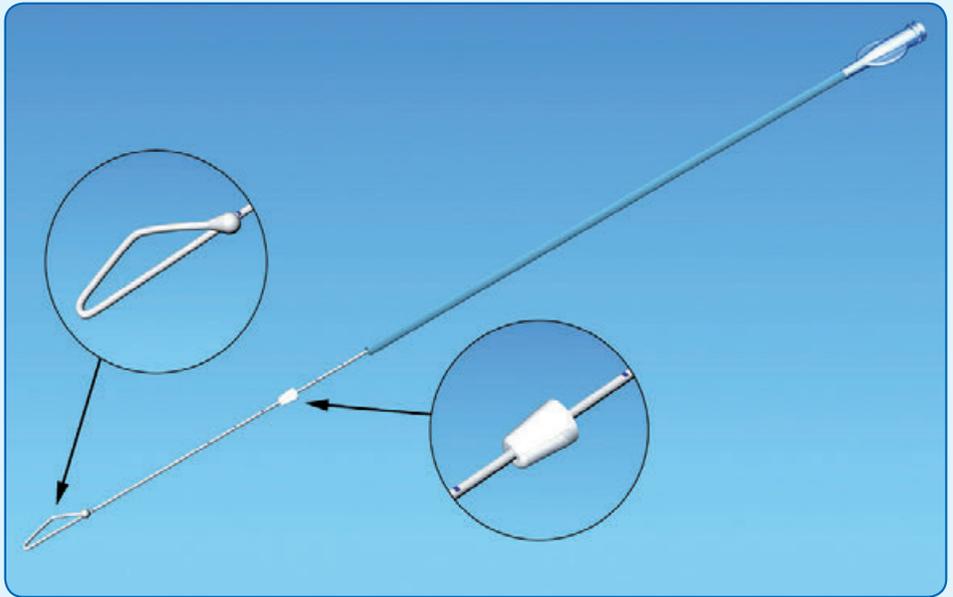


IUD EXTRACTOR

Dispositivo per la rimozione della spirale
Fixed IUD removal device

Cod.	Descrizione - Description
770123	lung. 28,5 cm - length 28.5 cm

Scatola da 10 pezzi. Fornito in confezione singola, monouso, sterilizzato a EtO e marcato CE.

10 pcs. per box. Supplied in single pouch, single use, EtO sterilized and CE marked.

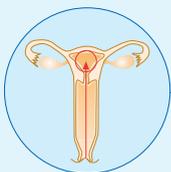

Indicato per la rimozione di IUD nei casi in cui non sia visibile il filo di repere o per le spirali che ne sono prive (tipo anello cinese). Il design della punta ne permette un facile inserimento e un agevole aggancio del filo o della struttura della spirale da rimuovere. L'olivetta al termine della punta evita traumi tissutali. La cannula in acciaio presenta una scala di misura di 9 cm (risoluzione 1 cm) ed è provvista di un fermo scorrevole per determinare la profondità di inserimento.

Designed for removing IUDs when the string is not visible or in presence of IUDs without strings (such as chinese ring). The design of the tip makes the instrument easy to insert and allows simplified hooking of the IUD string or structure. The ball at the end of the tip prevents tissue trauma. The steel cannula has a 9 cm graduation marking (resolution 1 cm) and is provided with a sliding positioner to determine the insertion depth.

