
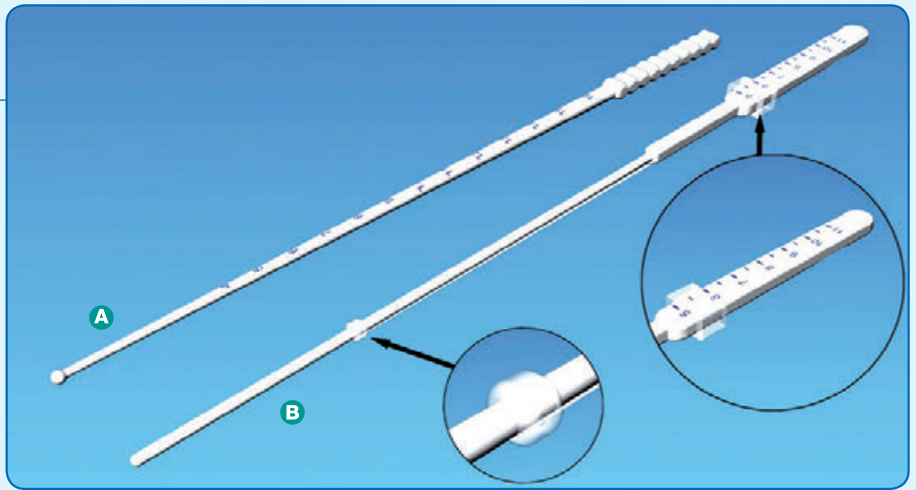


Endometrial sampling

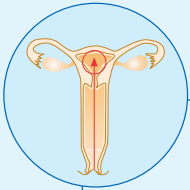

A RAMISTER

Isterometro retto
Straight Hysterometer

Misura la profondità dell'utero. Dritto, leggermente conico (da 3 a 5 mm di diametro) e graduato. La punta arrotondata permette un inserimento confortevole attraverso il canale cervicale e la malleabilità del materiale plastico consente di essere curvato manualmente a seconda delle necessità. La scala di misura va da 4 a 15 cm (risoluzione 1 cm) con lettura direttamente sulla sonda.

Used to measure uterus depth. It is straight, slightly conical (from 3 to 5 mm o.d.) and centimeter marked. Rounded tip allows a comfortable insertion through the cervical canal and malleable plastic permits manual bending of the instrument where necessary. Measurement scale ranges from 4 to 15 cm (1 cm resolution) and allows a direct reading on the probe.

Cod.	Descrizione - Description
720108	lung. 24,5 cm - length 24.5 cm

Scatola da 150 pezzi.
Fornito in confezione singola, monouso, sterilizzato a EtO e marcato CE.

150 pcs. per box.
Supplied in single pouch, single use, EtO sterilized and CE marked.

B RAMISTER-C

Isterometro curvo
Curved Hysterometer

Misura la profondità dell'utero. Ricurvo e leggermente conico (da 3 a 4 mm di diametro) per un inserimento confortevole attraverso il canale cervicale. Facile posizionamento grazie al fermo scorrevole da collocarsi sull'orifizio cervicale. La scala di misura va da 5 a 11 cm (risoluzione 0,5 cm) con lettura diretta sull'impugnatura.

Used to measure uterus depth. Curved shape and slightly conical configuration (from 3 to 4 mm o.d.) allow comfortable insertion through the cervical canal. It is easy to place thanks to the adjustable positioner to be set on the cervix external orifice. Measurement scale ranges from 5 to 11 cm (0.5 cm resolution) with direct reading on the handle.

Cod.	Descrizione - Description
720001	lung. 24,5 cm - length 24.5 cm

Scatola da 50 pezzi.
Fornito in confezione singola, monouso, sterilizzato a EtO e marcato CE.

50 pcs. per box.
Supplied in single pouch, single use, EtO sterilized and CE marked.

