


Catheters for Sonohysterography, Hysterosalpingography and Laparoscopy

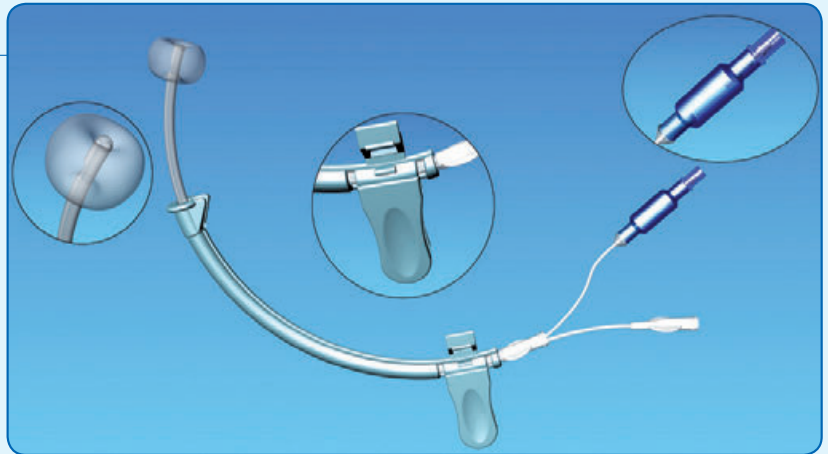

HANDLIN

Manipolatore/iniettore uterino
Uterine manipulator / Injector

- ▶ Impiegato in procedure dove viene richiesta la manipolazione dell'utero e/o per iniettare liquidi di contrasto (quali laparoscopia diagnostica e operativa, esami di fertilità, etc.). È caratterizzato da una forma ricurva che ne facilita l'inserimento. Provvisto di:
 - un catetere interno a due lumi con scala di misura da 4 a 10 cm (risoluzione 1 cm) per determinare la corretta profondità di inserimento;
 - un'impugnatura rigida;
 - un sistema di bloccaggio semplice e sicuro che impedisce movimenti accidentali durante la procedura;
 - un palloncino in PVC da 10 cc e 20 mm di diametro alloggiato sulla punta del catetere in modo da minimizzare la possibilità di perforazione accidentale ed evitare il reflusso del mezzo di contrasto;
 - una valvola unidirezionale per gonfiare il palloncino, assemblata a un polmoncino di riferimento che permette di visualizzare la reale condizione del palloncino interno;
 - un raccordo luer lock per l'iniezione del mezzo di contrasto.

- ✚ It is used for procedures where uterus manipulation and/or injection of contrast media is requested (e.g. diagnostic and operative laparoscopy, fertility examinations, etc.). It has a curved shape for an easy insertion. It is composed by:
 - double lumen inner catheter with graduation from 4 to 10 cm (1 cm resolution) to assess correct insertion depth,
 - stiff handle allowing a firm and safe manipulation;
 - simple locking system to avoid accidental movements during the procedure;
 - 10 cc PVC balloon with diameter 20 mm placed on the catheter tip to minimize any accidental perforation and avoid contrast media reflow;
 - one-way inflation valve directly connected to a pilot balloon showing the real condition of the inner balloon;
 - luer lock connector for contrast media injection.

Cod.	Descrizione - Description
720199	Ø esterno 5 mm - O.D. 5 mm

Scatola da 10 pezzi.
Fornito in confezione singola, monouso, sterilizzato a EtO e marcato CE.

10 pcs. per box.
Products are single packed, single use, EtO sterilized and CE marked.

